

**From:** EMIM 2017 abstract@eventclass.com  
**Subject:** Abstract transmission report I EMIM late-breaking abstract 2017  
**Date:** 22 February 2017 at 23:41  
**To:** Greetje Vande Velde greetje.vandevelde@kuleuven.be  
**Cc:** Emiliano D'agostino emiliano.dagostino@dosevue.com, Christian Vanhove Christian.Vanhove@UGent.be

E

**Dear Prof. Dr. Vande Velde**

You have successfully submitted your abstract for the EMIM 2017. This is the transmission report. You can make changes on your submission until the deadline. To modify your submission you need to log in to the control center. The login data for the control center has been sent to the corresponding author in a separate e-mail.

**Your registration:** With the submission of a late-breaking abstract you are entitled to register for the early bird fee. Deadline for early bird registration was February 14, 2017. In case you have not registered yet, please enter your individual reduction code **d2017JUSze** during your online registration. This code is valid until March 27, 2017. Registration is possible via [www.emim.eu](http://www.emim.eu). In case you are not the presenting author of this submission please forward the code to the presenting author.

Looking forward to meeting you in Cologne. For further information on programme, speaker, travel... [www.emim.eu](http://www.emim.eu).

**Data:**

ID: **#516**  
Topic: **Imaging Technology, Methodology & Data Analysis > Nuclear and CT Imaging**  
Type of presentation: **Oral or Poster presentation**  
Keywords: **dosimetry, micro-CT, radiotherapy**

**Authors:**

Main author: **D'agostino Emiliano**  
Presenting author: **D'agostino Emiliano**  
Correspondence author: **Vande Velde Greetje**  
Authors including in correspondence: **D'agostino Emiliano, Vanhove Christian**  
Access rights\*: **Vande Velde Greetje**  
**<[greetje.vandevelde@kuleuven.be](mailto:greetje.vandevelde@kuleuven.be)>**

\*This author already has access rights for control center or the access data were sent automatically by the system to this author.

**Your submission:**

## **DoseWire: optical fiber-based real-time dosimetry for microCT and radiation therapy applications**

**D'agostino Emiliano**<sup>1</sup>, Bogaerts Ria<sup>2</sup>, Vanhove Christian<sup>3</sup>, Vande Velde Greetje<sup>4</sup>

<sup>1</sup>DoseVue NV - Mol, Belgium

<sup>2</sup>Laboratory of Experimental Radiotherapy, Oncology - Leuven, Belgium

<sup>3</sup>Institute Biomedical Technology, Medical Imaging and Signal Processing - Gent, Belgium

<sup>4</sup>Biomedical MRI/ MoSAIC, Imaging and Pathology - Leuven, Belgium

### **Introduction**

Ionizing radiation is frequently used in pre-clinical research for diagnostic and therapeutic purposes. Radiation based imaging, such as micro-computed tomography ( $\mu$ CT) is a powerful tool for non-invasive longitudinal follow-up of animal models. A lot of interest has also been directed towards the optimization of small animal radiation therapy (SmART), for subsequent translation to patients. In both cases (imaging and therapy), the animal models involved in the studies may be exposed to significant amounts of radiation. Therefore, an accurate monitoring of radiation doses is crucial to avoid excessive exposure of healthy tissues.

DoseWire is a dosimetry device based on an optical fiber, allowing real-time radiation monitoring. The device is portable and user-friendly and can also be used within a mouse for in-situ measurements in a minimally invasive way. DoseWire is characterized by an elevated signal to noise ratio and can therefore be used without need of corrections for the so-called stem effect [1]. This is particularly important when using large irradiation fields (as in imaging) and/or when the beam axis is close to 45°.

## Methods

DoseWire was used with a small animal  $\mu$ CT (SkyScan 1076, 1276 and 1278, Bruker microCT, Kontich, Belgium) and on a small animal irradiation platform (SARRP, Xstrahl, Surray, UK). Measurements on the  $\mu$ CT were performed for different protocols (changing mAs, kV (35-70 kV range) and filter combinations) and scanner types. Measurements on the SmART system were all performed at 220keV using different tube currents (dose rates). In all cases, DoseWire was calibrated by an ionization chamber.

## Results

The graph below shows a linearity test for measurements performed on a SmART system, demonstrating that the response of the DoseWire is linear within a broad range of dose rates. Measurements were reproducible with a variability below 1-2%. DoseWire was also tested on imaging fields, using different protocols covering a broad range of microCT imaging applications on different scanner types. DoseWire allowed real-time feedback on the effect on dose rate when changing scan parameters (tube current, filter type, kV), thereby allowing optimization of scan settings in function of minimal dose. Reproducibility of the measured dose within the FOV, to the animal and to organs was found to be excellent.

## Conclusions

Real-time dose measurements are an important asset when performing studies involving intensive use of ionizing radiation, either for imaging purposes, or for therapy studies. In both cases, monitoring the dose in real-time can contribute to optimized experimental protocols, thereby allowing to maximize the quality of the study and to prevent over-exposure of tissues, possibly leading to significant experimental bias. DoseWire proved to be an accurate and user-friendly tool for systematic real-time dose monitoring in different experimental conditions.

[1] On the nature of the light produced within PMMA optical light guides in scintillation fiber-optic dosimetry, Theriault-Proulx et al. *Phys. Med. Biol.* 58 2073

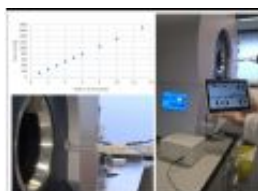


Figure 1:

**Figure 1.**  
DoseWire

Best regards

Your EMIM 2017 abstract management team

---

**Scientific contact:**

ESMI - European Society for Molecular Imaging  
Management Office  
Doris Kracht  
phone: +49 221 78966180  
e-mail: [office@e-smi.eu](mailto:office@e-smi.eu)

**Technical support:**

EventClass GmbH  
phone: +49 351 30900031  
e-mail: [abstract@eventclass.com](mailto:abstract@eventclass.com)